


TSIII Ø3.0 Implant

Strong Mini Solution

OSSTEM[®]
IMPLANT




Implant designed for safe use in extremely narrow bone, where previous two-piece implants are difficult to install

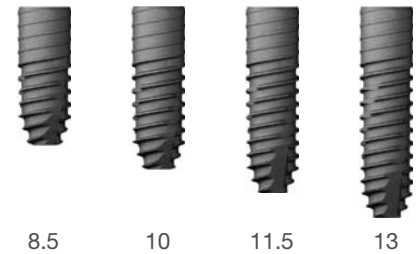

Smaller but stronger fatigue strength on fatigue life

Fatigue life (cycle)

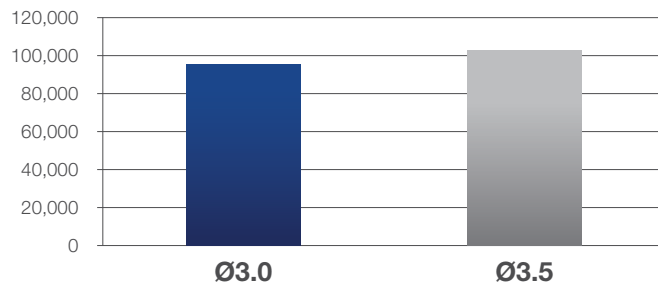


Fracture risk is reduced by reinforcing fatigue strength, to a level similar to that of Ø3.5 fixture

※ 5mm cantilever fatigue test (applied loading : 350Ncm)

100% compatible with upper TS mini



- 100 % compatible with upper TS mini
- Requires the cover screw and analog for Ø3.0 only

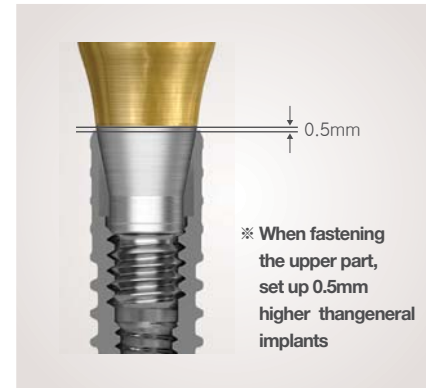


Drilling Sequence & Spec.

Drilling Sequence



Bone quality	Ø 2.2 drill	Ø 2.7 drill	F3.0 cortical drill	3.0 fixture
Soft	●			
Normal	●	●		Implant placement
Hard	●	●	●	

Spec.

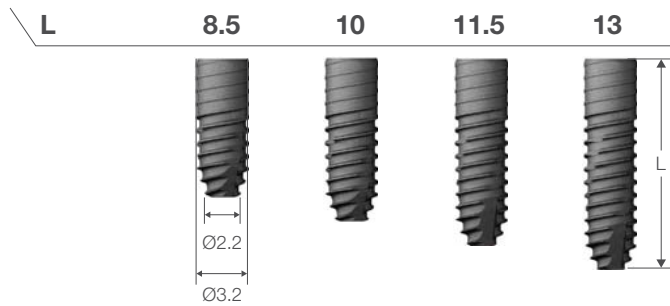

OSSTEM[®]
IMPLANT

DEUTSCHE OSSTEM GmbH

Mergenthalerallee 25, 65760 Eschborn, Germany
Tel: +49(0)6196 777 550 Fax: +49(0)6196 777 5529
E-mail: info@osstem.de www.osstem.de